LU 8.6 Twists in the fibres

Introduction

Has anybody in your family bought silk or wool and found later that it was not pure silk or wool? Although we use bags, ropes, clothes, and items made from different fabrics, it is difficult to ascertain the purity of the fabric by touch, texture, or weight alone. A fibre or yarn may look like cotton but it may be synthetic. One cannot trust appearance alone.

Would you like to learn simple ways to identify fibres? In this Learning Unit, you can explore this through two techniques: burning test and microscopy.

Did you know?

A thread is not a fibre but a bundle of fibres. You can see this in figure 1.

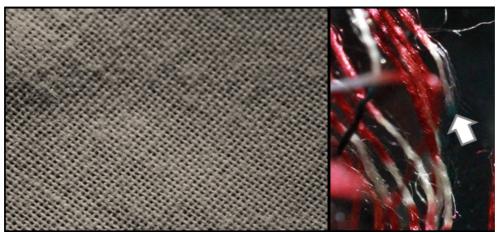


Figure 1: (Left) A cloth is made of threads. (Right) Each thread is made of several fibres (indicated by arrow).

Materials required

- Threads: cotton, polyester (or any other synthetic thread), wool (if you don't get real wool, a broken hair can also be used), and a thread made of unknown fibre. The fibres should be clean, so that any other substance sticking to them does not affect the observations.
- Candle, match stick, beaker, forceps, watch glass, water, tray.
- Microscope, glass slide, and cover slip.

Burning test for fibres

There is a popular saying in Hindi "Rassi jal gai, par bal nahin gaye.", meaning "the rope got burnt, but the twists in the rope have remained (in the form of ash)". This is a property of ropes that are made of plant fibres, which on burning produce ash in which twists can still be seen. Figuratively, it also refers to someone's personality traits that did not change even after facing a lot of difficulties in life.

Textile experts across the world have been using burning tests to identify fibres. Plant fibres burn slowly producing ash, which glows for few seconds (known as afterglow) before cooling down. The ash roughly retains the shape of fibres but can be easily crushed to powder.

Animal-based fibres such as silk and wool also burn to give ash but in much lesser quantity. Sometimes they burn producing a small bead which gets crushed easily. Moreover they burn with

an odour, similar to that of burning meat.

Synthetic fibres usually burn much faster than plant or animal fibres. When exposed to flame, they melt and then decompose, producing a lot of heat. These do not form any ash but form a bead which is not easy to crush.

Safety Precautions: Be very careful while doing the burning tests. Do not be too close to the flame or the burning fibre and do not throw burnt threads around. Some fibres burn very vigorously and their melts can also cause burns or damage to surrounding objects. Keep water in a watch glass or a container nearby, and put the burnt fibres in the water.

Fibres under the microscope

Different fibres have different shapes and surface features that can be seen under a microscope. Depending on the conditions in which the fibre has formed (shape of original plant cell or animal cell(s) and how the fibre was dried), it can have a circular, elongated, or irregular cross section. Cotton fibres, in particular, have a non-circular cross-section due to which the twists in the fibres can be easily observed under a microscope. Plant fibres generally have rough surfaces. In a bundle of natural fibres, fibre thickness may also vary from one fibre to another.

Synthetic fibres usually have circular cross-sections and smooth surfaces because they are formed by passing molten polymer through circular holes, similar to the way noodles are made. Thickness of a synthetic fibre is uniform along the length of a fibre, and is also the same for different fibres in a thread.

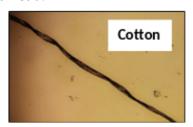






Figure 2: Fibres under microscope

Animal fibres such as wool and silk have circular cross-sections and surfaces smoother than plant fibres but rougher than synthetic fibres. Wool/hair also have scales of the surface (which may not be observed if the wool/hair has been treated with strong bleaches, dyes, or chemicals that damage the hair scales).

However, all fibres on use (even synthetic fibres) usually develop surface roughness due to wear and tear.

Q1. What differences do you observe between the microscopic images of cotton and polyeste
fibres in figure 2? Are there any other different features of these fibres that you can conclude from
these observations?

Task 1: Known fibres

Take a cotton and a polyester or synthetic thread (like nylon or acrylic).

1. Observe the fibres in each thread with the naked eye. Note their physical properties such as shiny, dull appearance, and if they are smooth or rough to touch.

Fibre	Observations (Shiny/dull, rough/smooth texture)

Table 1

2. Take a tray, and fix a candle in the middle of it. Fill the tray slightly with water. Light the candle. Hold the thread with a pair of forceps or tongs and bring one end of it close to the flame. Note your observations about the thread burning in the following table. Collect the ash/bead formed on a watch glass. [Caution: Keep your head/body parts away from flame as some fibres burn very vigorously. Extinguish flames of any burning fibres in water; do not throw them anywhere else.]

Note the following observations for each fibre.

Sr. No.	Did it melt? (Yes/No)	Any smoke? (Yes/No)	Smell (like burning paper or plastic)	Ash/ Bead formed	Afterglow? (Yes/No)

Table 2

3. Take some water in a glass. Check if it is acidic or basic with a litmus paper. Add one drop of water to the ash/bead on the watch glass. Wait for 1-2 minutes and check with litmus papers (red and blue)/a drop of phenolphthalein solution/a pinch of turmeric. Has the water become acidic or basic on contact with the ash/bead? Note: Ash or bead will not completely dissolve in water; some solid will remain in both cases.

Sr. No.	Thread burnt	Ash + Water (Neutral/Basic/Acidic)

Table 3

4. Now observe the thread under a microscope as described below. (Did you know that Antonie van Leeuwenhoek, considered to be a pioneer of microbiology, started his journey as a colth merchant. He started using self-made lens arrangements to observe the quality of fibres in his shop. This setup of lenses was a rudimentary microscope which later also helped him to study animals and organisms in detail).

A thread usually has several fibres bundled together. Using a pin or forceps, loosen out the fibres in a thread and pull out a fibre. Put the fibre on a glass slide and cover it with a cover slip. Observe it under 10X objective. If you are not able to focus on the fibre, make sure that the fibre is under the objective lens and while focusing, keep the distance between the lens tip and the cover-slip is around 0.5 cm.

Note the features of the fibre. Next put 3-4 fibres together on the slide and observe the variation in the thickness of different fibres.


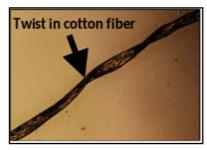


Figure 3

Q2. Is the fibre that you observe uniform in thickness along the length or is the thickness different at middle and ends of the fibre (To see this, you will to move the slide on the stage to see it from one end to another end)?

- Q3. Do you observe twists in the cotton fibres?
- Q4. Do you observe folds in the cotton fibre?
- Q5. What can you say about the thickness of various fibres?

Now sketch the shape of the fibre observed under microscope, and write its features:

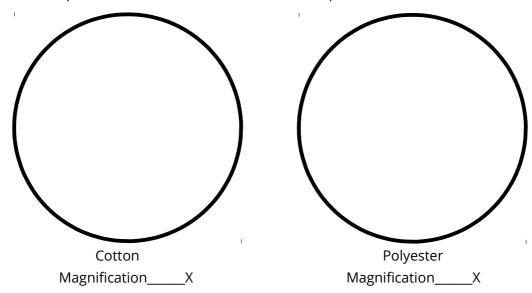

Fibre	, ,	Are all the fibres same or different in shape?
Cotton		
Polyester		

See the sketches done by your classmates and copy here at least one different sketch of cotton and polyester as drawn by one of your classmates .

Cotton:	Polyster:
Q6. Based on your observations, what can you conclud plant-based fibre), and polyester (a synthetic fibre).?	e about the features/properties of cotton (a
Q7. Why do you think the ash obtained from burning water, and the bead from polyester did not?	cotton changes the (acidic/basic) nature of


Task 2: Wool/hair (animal fibre)


Take a woollen thread or human hair (because real wool is also the hair of some animal) and perform the following steps.

1. By burning test as done in Task 1, check if it is natural or synthetic. (If it is synthetic, then try finding a real wool/hair sample.)

How did it melt (smell, smoke, ash/bead formed)?	afterglow,		
Ash/bead + water: Acidic or basic or neutral?		Crushable?	

2. Separate fibres out of the thread, as done in Task 1 and observe the fibres under microscope.

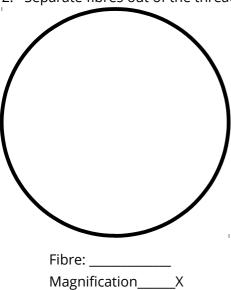

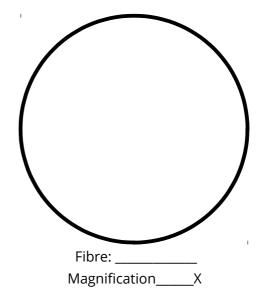

ibre description (straight or twisted, transparent or opaque, uniform or variable thickness along enght)
are all fibrac came or different in chang?
Are all fibres same or different in shape?

The science behind the burning tests

Now that you have done some burning tests, let us try to understand why plant-based, animal-based, and synthetic fibres burn differently.

- Plant-based fibres (cotton, linen, hemp, sugarcane, jute) have high amount of cellulose which burns easily. However, these also have some amount of silica and metals (~0.1-0.6% by mass) which lead to ash formation (~0.5-1.2% by mass). These fibres burn slowly with a flame. When the flame goes off, the remaining silica and metal compounds glow red like burning coal. As the burnt fibre/ash cools down, the afterglow disappears and fine ash is obtained.
- Synthetic fibres usually do not have silica or metals. These are often produced from polymers made of non-metallic elements such as carbon, nitrogen, oxygen, and hydrogen. These fibres have low melting temperatures. Therefore when exposed to flame, these melt and then decompose. Burning a synthetic fibre produces a lot of heat but may or may not form smoke (Smoke formation depends on percentage of carbon in the polymer. If carbon percentage is high, then all of carbon is not able to get enough oxygen to form carbon dioxide, and hence particles of unburnt carbon and related substances are formed which become smoke).
- Animal-based fibres (wool, silk) are predominantly made of proteins (such as keratins) which burn with a smell, similar to obtain from burning of other proteins such as meat.

Task 3: Unknown fibre

1. Take a thread of unknown material. By observing this thread with the naked eye, guess if it is a natural or a synthetic fibre.

Sr. No.	Observations	Natural/Synthetic

Table 5

- 2. Conduct the burning test for the thread as done in Task 1, identify if it is natural or synthetic. [Remember the precautions of keeping your head/body parts away from flame as some fibres burn very vigorously. Extinguish flames of any burning fibres in water; do not throw them anywhere else.]
- 3. Separate the fibres out of the thread, as done in Task 1, and observe the fibres under a microscope.

Sketch	Describe a single fibre	Variation in a bunch of fibres

Table 6

Note: Do check if the thread you have is not mixed (i.e., if it consists of more than one kind	of
fibres) and hence may give properties of both kinds of fibre, such as giving both ash and bead	on
burning.	

Based	on	the	above	tests,	try	to	identify	the	fibre:	The	observed	l fibre	is
			beca	use									
4. Paste	e/stap	ole a s	ample of	the thre	ead he	ere.		1					
As a cu liveliho Q1. For	ustom ods o each	ner, the f million of the	ons of pe e fibres t	e of fibr ople de hat you	pend identi	on tl fied,	he fibres t , list the pe	nat we eople v	choose who are	to use involve	nought abo ? ed from the hen you bu	product	
	or any		-				-		-		sed for makepare the th	_	
Q3. Bas	sed o	n wha	t you hav	e learne	ed abo	ove,	can you sa	ay if th	e wicks	for can	idles and o	il lamps o	can

be made using synthetic fibres? Why?

Vigyan Pratibha Learning Unit	Twists in the fibres			
Q4. Use in wick-making continues to be an important reason fo professions that depend on the use of cotton in wicks.	r the sale of cotton. Name the			

Credits

Main Author: Ankush Gupta

Contributing Authors: Chitra Joshi, Bhagyashree Thube, Krupa Subramaniam

Reviewers: Savita Ladage, S. D. Samant, Madanrao D.

Editors: Beena Choksi, Geetanjali Date, Ankush Gupta, Reema Mani, K. Subramaniam

Creative Commons Licence: CC BY-SA 4.0 International, HBCSE